



Hearing Systems – Part 2 of 2

<https://hearingsystems.github.io/>



TU Ilmenau – Audio Systems Technology

(Online) | 27. Januar 2021

Dr.-Ing. Tamas Harczos


•• Reminder

- Useful links:
 - <https://hearingsystems.github.io/>,
 - <https://www.tu-ilmenau.de/mt/lehrveranstaltungen/lehre-fuer-master-mt/audio-systems-technology/>.
- Books:
 - *Hörakustik 3.0 - Theorie und Praxis* (Jens Ulrich, Eckhard Hoffmann), ISBN: 978-3-9428-7336-9
 - *Auditory Prostheses: New Horizons* (Zeng, Popper, Fay, eds.), ISBN: 978-1-4419-9434-9
 - *Human and Machine Hearing* (Richard F. Lyon), ISBN: 978-1-1070-0753-6, see [draft online](#)
- Programming / Testing:
 - [Oldenburg openMHA](#) (Master Hearing Aid), [Hearing Loss Simulator](#) (CDC, USA)
 - [Cochlear implant \(CI\) and hearing loss simulator](#) (AngelSim)
 - [CI simulation](#) (UT Dallas, USA), [CI simulation](#) (UGR, Spain)
 - [MATLAB](#) 2018b+ (free student version) or [GNU Octave](#) 6.1+ (free)
 - [Wavosaur](#) or [Audacity](#) audio editor (both free)

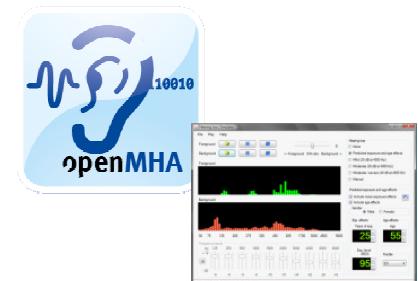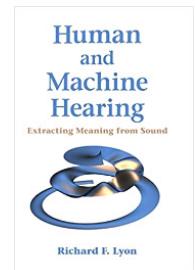

•• Reminder

- Useful links:

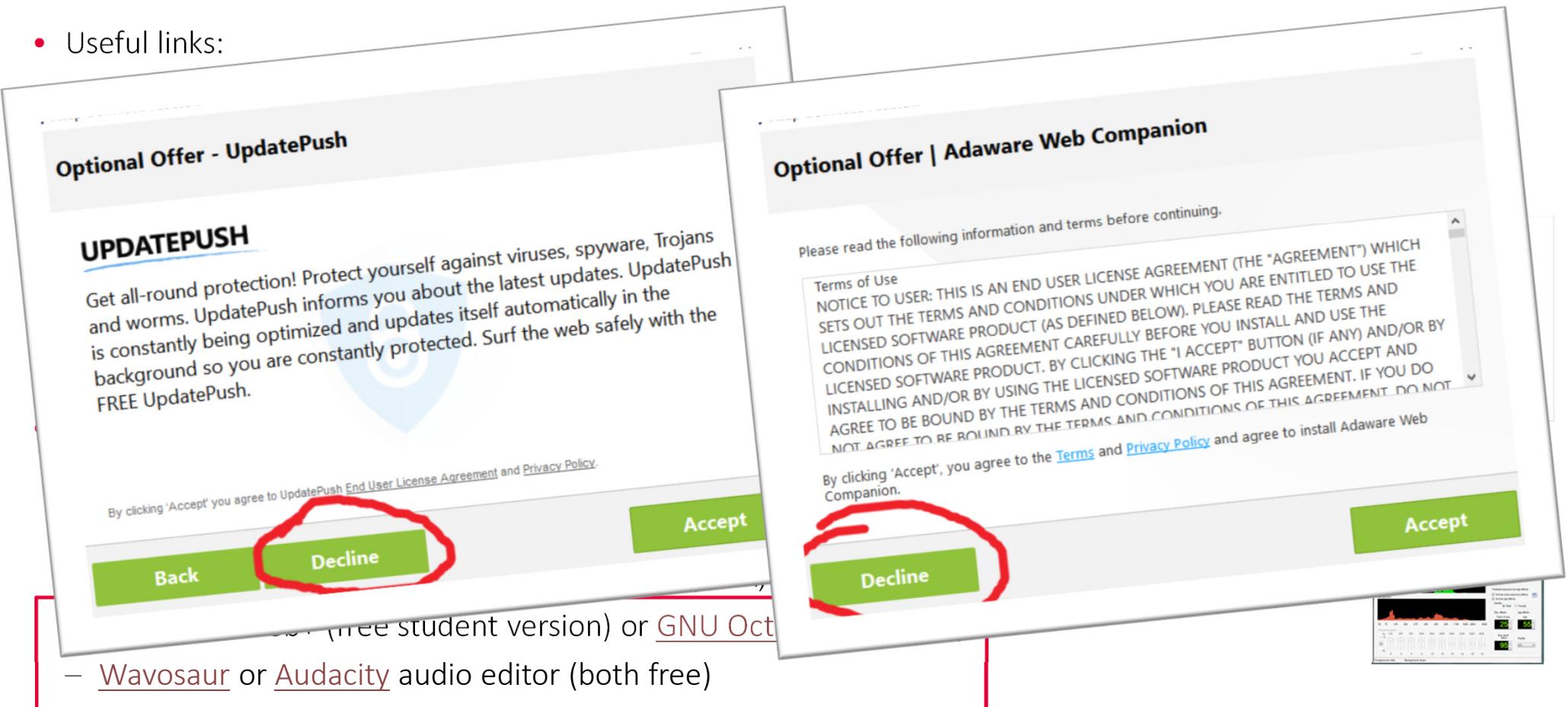

Optional Offer - UpdatePush

UPDATEPUSH

Get all-round protection! Protect yourself against viruses, spyware, Trojans and worms. UpdatePush informs you about the latest updates. UpdatePush is constantly being optimized and updates itself automatically in the background so you are constantly protected. Surf the web safely with the FREE UpdatePush.

By clicking 'Accept' you agree to UpdatePush End User License Agreement and [Privacy Policy](#).

Accept **Decline** (circled in red)

– [Audacity](#) (free student version) or [GNU Octave](#) (free) – [Wavosaur](#) or [Audacity](#) audio editor (both free)

Optional Offer | Adaware Web Companion

Please read the following information and terms before continuing.

Terms of Use
NOTICE TO USER: THIS IS AN END USER LICENSE AGREEMENT (THE "AGREEMENT") WHICH SETS OUT THE TERMS AND CONDITIONS UNDER WHICH YOU ARE ENTITLED TO USE THE LICENSED SOFTWARE PRODUCT (AS DEFINED BELOW). PLEASE READ THE TERMS AND CONDITIONS OF THIS AGREEMENT CAREFULLY BEFORE YOU INSTALL AND USE THE LICENSED SOFTWARE PRODUCT. BY CLICKING THE "I ACCEPT" BUTTON (IF ANY) AND/OR BY INSTALLING AND/OR BY USING THE LICENSED SOFTWARE PRODUCT YOU ACCEPT AND AGREE TO BE BOUND BY THE TERMS AND CONDITIONS OF THIS AGREEMENT. IF YOU DO NOT AGREE TO BE BOUND BY THE TERMS AND CONDITIONS OF THIS AGREEMENT, DO NOT

By clicking 'Accept', you agree to the [Terms](#) and [Privacy Policy](#) and agree to install Adaware Web Companion.

Accept **Decline** (circled in red)

•• Recap: Hearing aids – Modern form factors



Behind the ear (BTE)



Receiver in canal (RIC)



In the canal (ITC)



Completely in canal (CIC)

Source: betterhearing.org, 2021



- + least earwax impact
- definitely visible



- + least feedback



- + glasses-friendly
- ear occlusion



- + glasses-friendly
- ear occlusion
- usability

•• Recap: Hearing aids – Features

General:


- Small, lightweight, intuitive, discrete.
- Remote control, smartphone app.
- Wireless streaming (TV, remote microphone).
- Quickly rechargeable, long lasting battery.

Audiological features:


- Low noise, high bandwidth (10 kHz+).
- 16+ frequency channels for DRC.
- Adaptive microphone directionality.
- Noise reduction: ambient, wind, impulse.
- Environmental classification.
- Advanced feedback cancellation.
- Frequency shifting.
- Low latency (<< 10 ms), binaural algorithms.

Fitness and medical features:


- Fitness tracking (step counter, pulse rate etc.).
- Blood pressure monitoring.
- On-demand activity log.
- Fall detection.

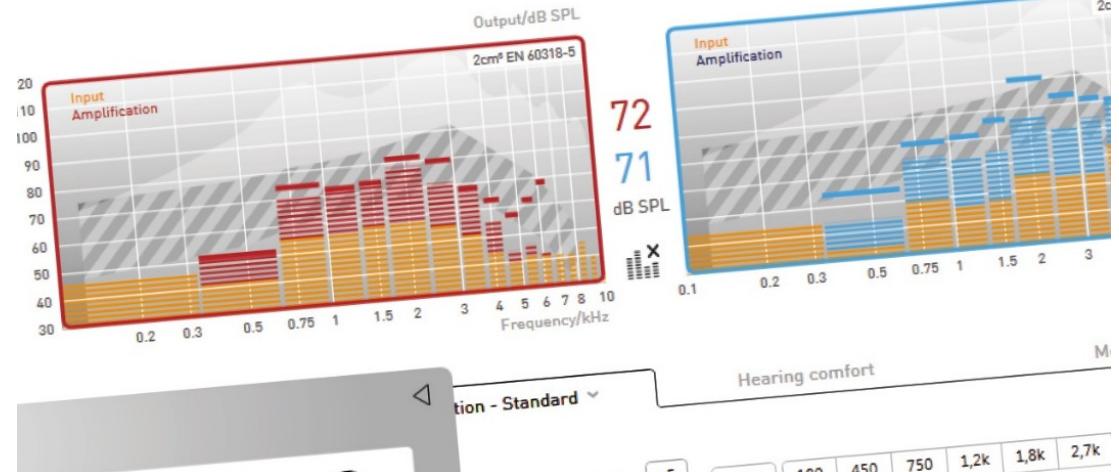

•• Hearing aids – What is inside?

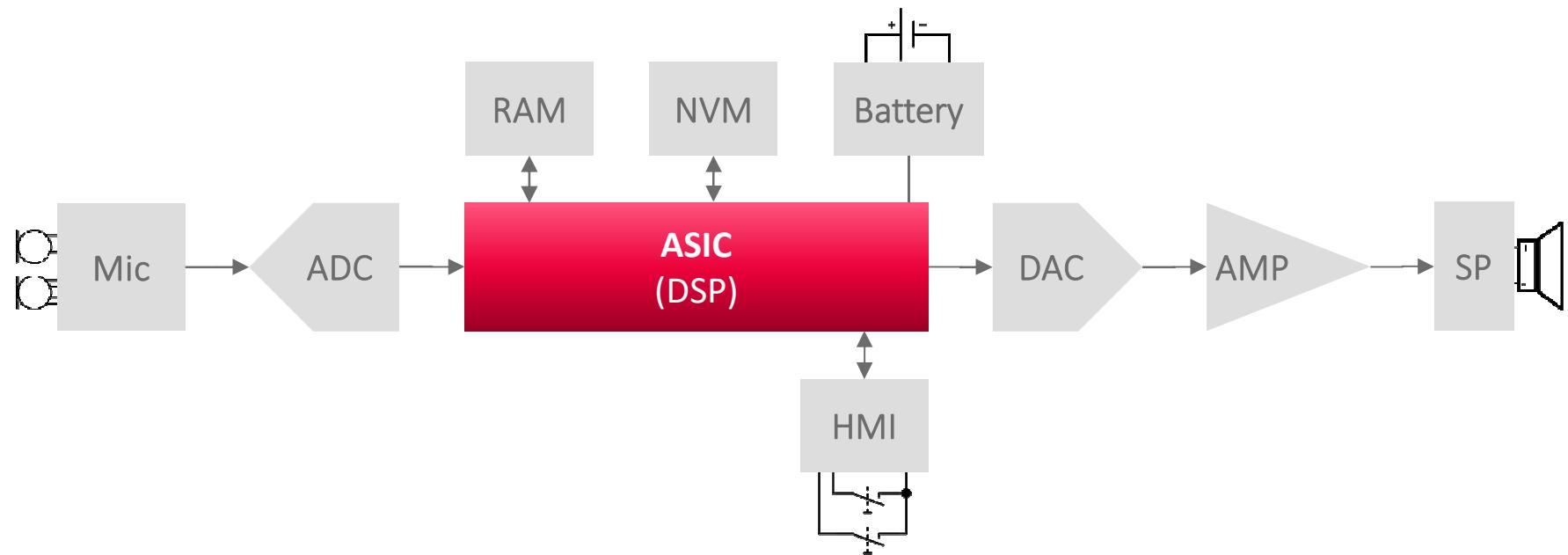


•• Hearing aids – What is inside?

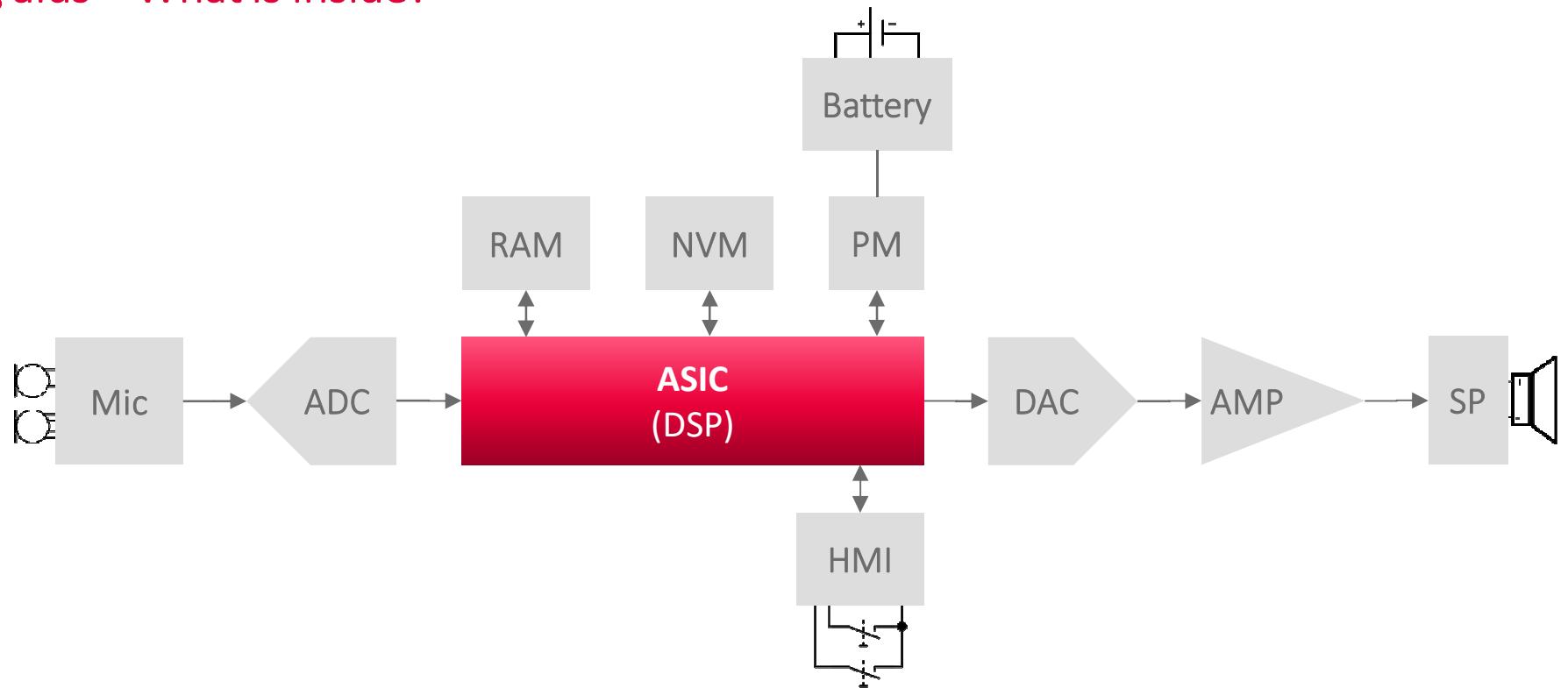

ASIC: application specific integrated circuit, **HMI:** human-machine interaction, **NVM:** non-volatile memory, **PM:** power management


page 7

•• Hearing aids – What is inside?

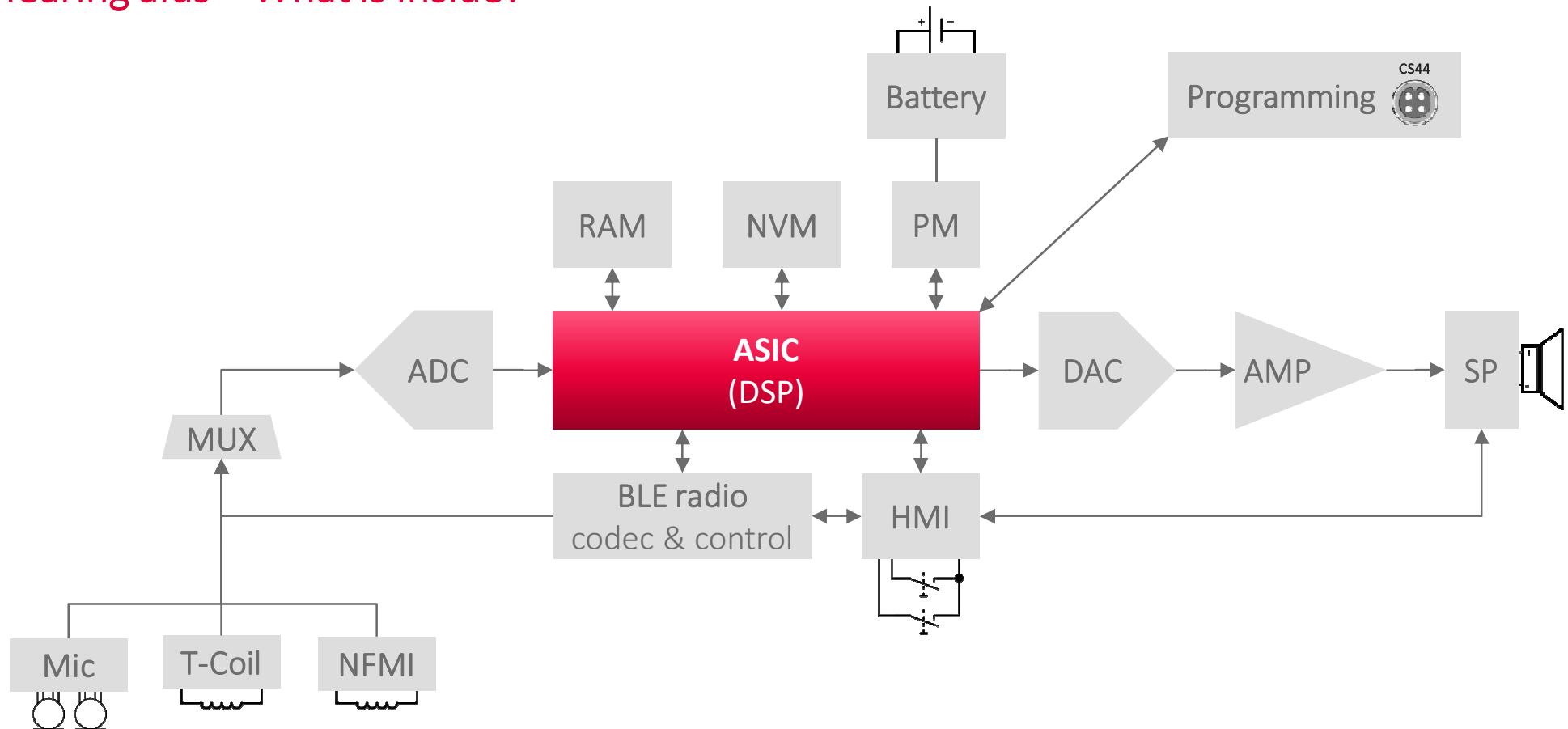


ASIC: application specific integrated circuit, **HMI:** human-machine interaction, **MUX:** multiplexer, **NFMI:** near-field magnetic induction, **NVM:** non-volatile memory, **PM:** power management, **T-Coil:** Telecoil (receiver for inductive loop), **WL:** (other) wireless technology

•• Hearing aids – What is inside?

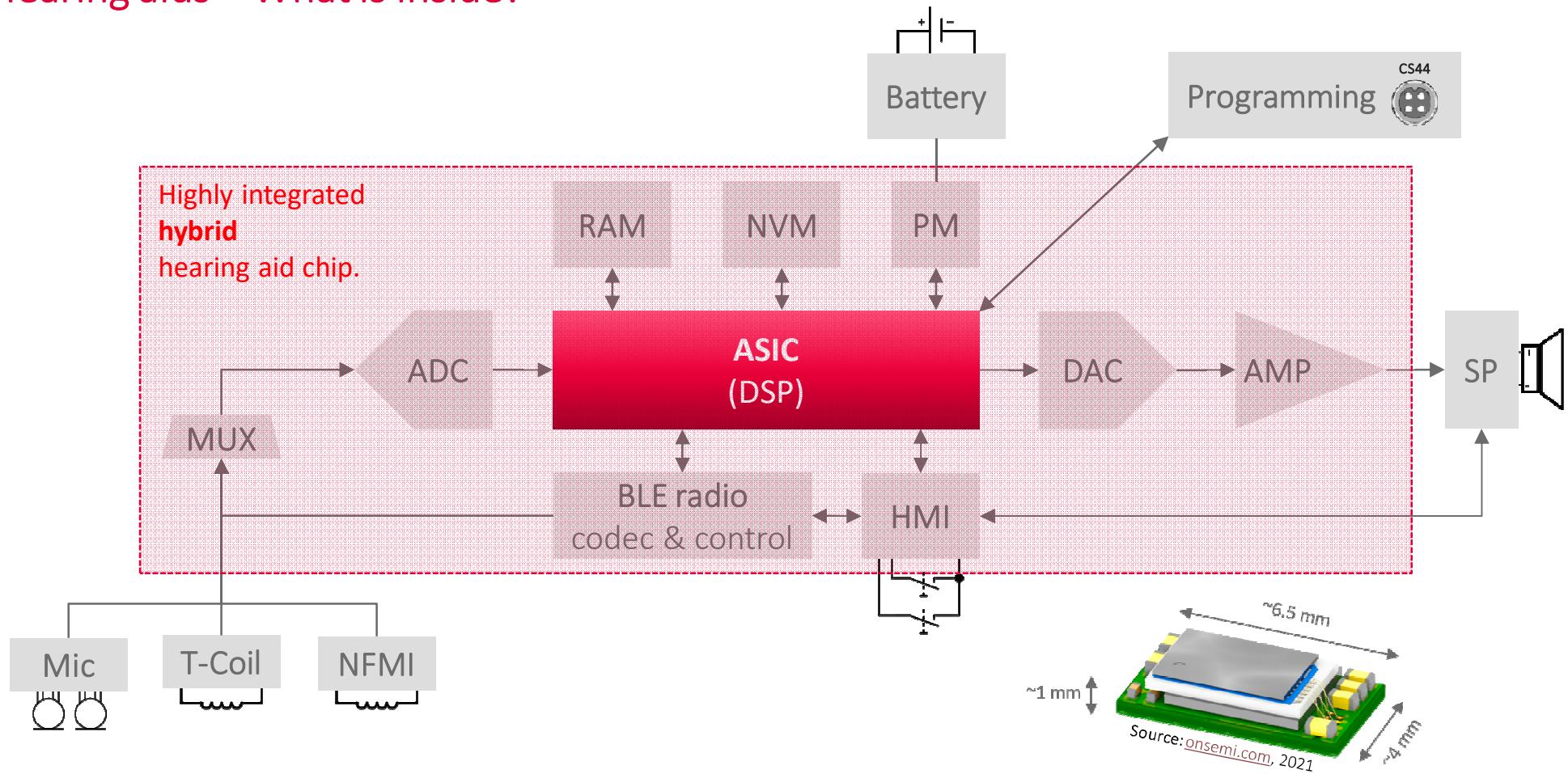


ASIC: application specific integrated circuit, **HMI:** human-machine interaction, **MUX:** multiplexer, **NFMI:** near-field magnetic induction, **NVM:** non-volatile memory, **PM:** power management, **T-Coil:** Telecoil (receiver for inductive loop), **WL:** (other) wireless technology

•• Hearing aids – Basic signal flow

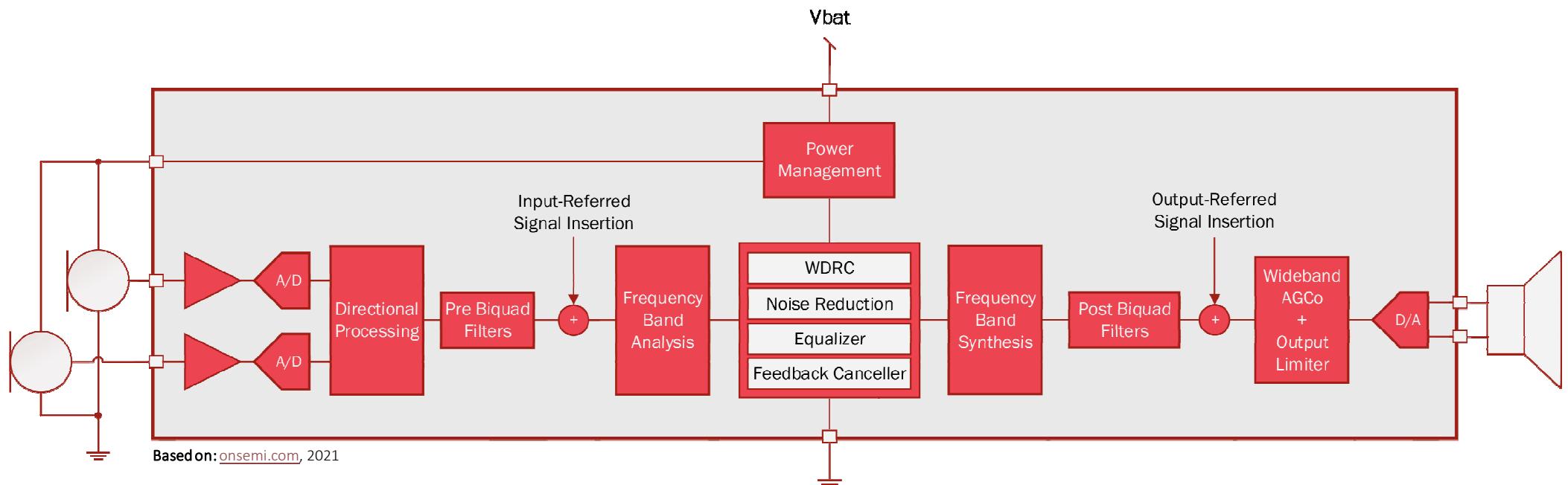

•• Hearing aids – Architecture of an ASIC

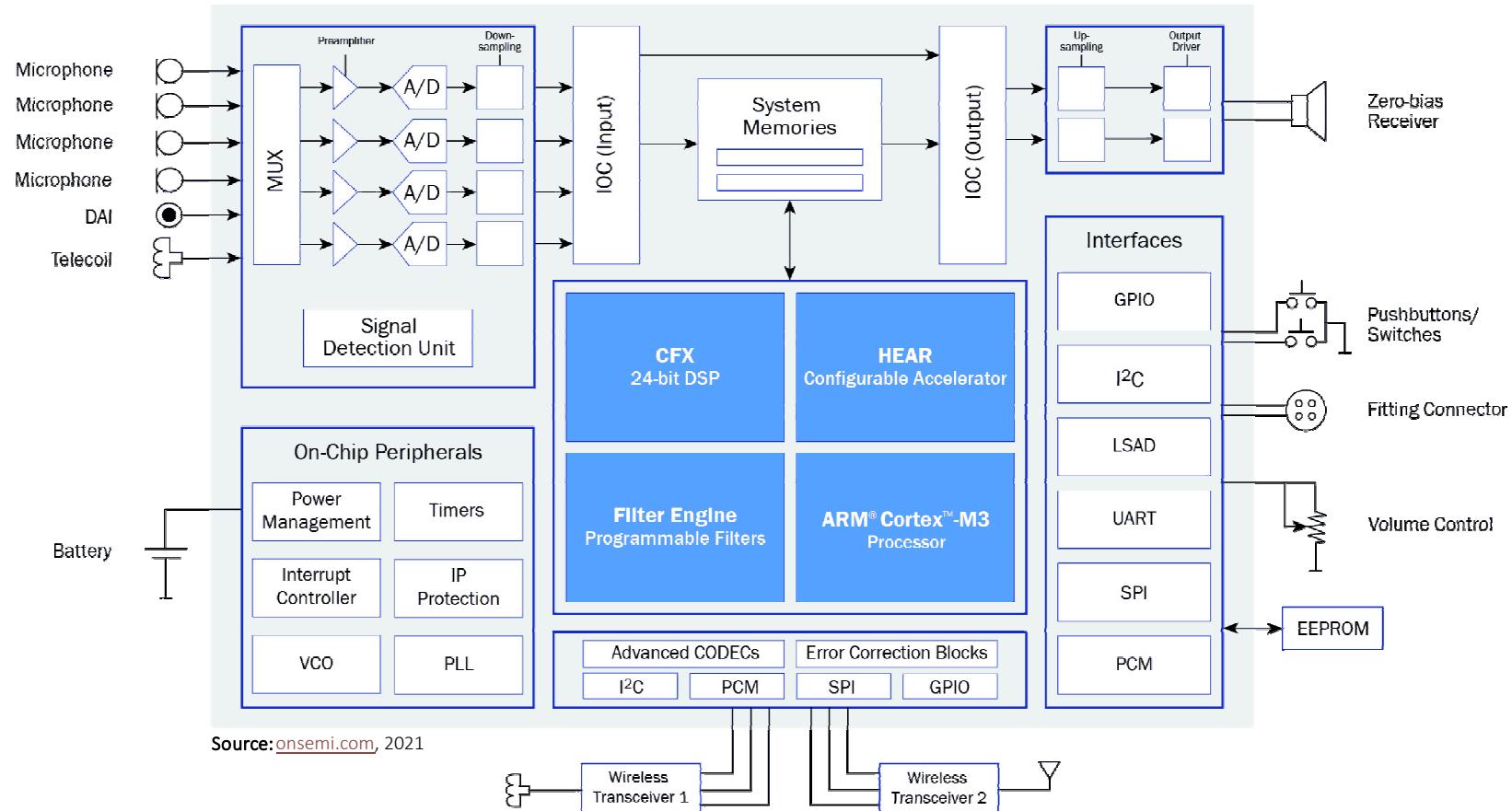

DAI: digital audio input, EEPROM: electrically erasable programmable read-only memory, GPIO: general purpose input/output, I²C: inter-integrated circuit, IOC: input/output controller, IP: intellectual property, LSAD: low-speed A/D converter, PCM: pulse code modulation, SPI: serial peripheral interface, UART: universal asynchronous receiver transmitter, WL: (other) wireless technology

•• Hearing aids – Basic specifications

Rough estimates!




Source: jabra.com, 2021

Source: amd.com, 2021


	Hearing aid ASIC	Hearable ASIC	PC CPU
Typical power consumption (algorithms on, wireless off)	~ 1-2 mW	~ 15-50 mW	~ 35-125 W
Typical power consumption (algorithms on, BLE in use)	+ 3-8 mW	+ 15-50 mW	
Clock speed	~ 10-60 MHz	~ 20-200 MHz	~ 3000-4000 MHz
MIPS @ max. clock speed	~ 500 MIPS	~ 2000 MIPS	~ 500 000 MIPS
Typical dimensions	~ 4 x 7 x 1 mm = 28 mm ³	~ 5 x 8 x 1 mm = 40 mm ³	~ 40 x 40 x 6 mm = 9.6 cm ³
Number of cores	~ 2-8	~ 4-8	~ 2-32
Fabrication process	~ 45-22 nm	~ 45-22 nm	~ 14-7 nm

•• Hearing aids – Construction of a hearing aid

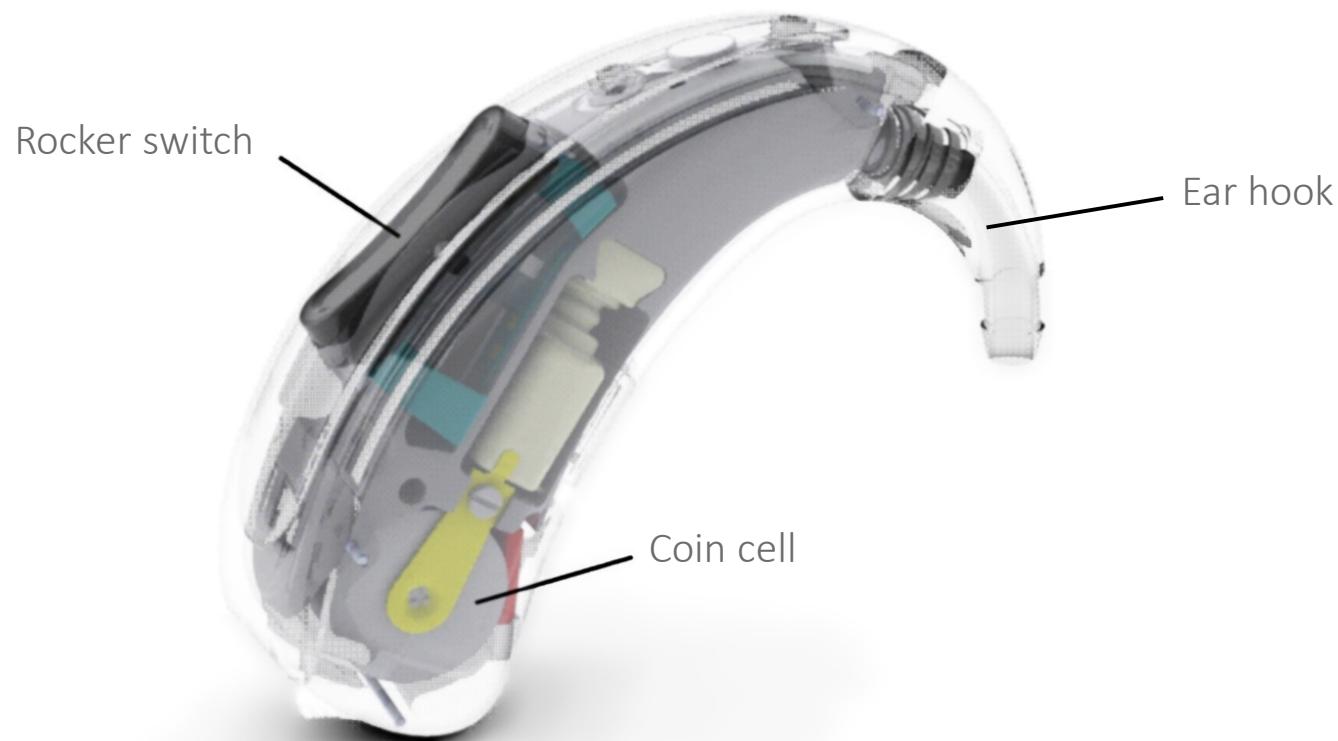




Source: [kind.com](#), 2021

•• Hearing aids – Construction of a hearing aid

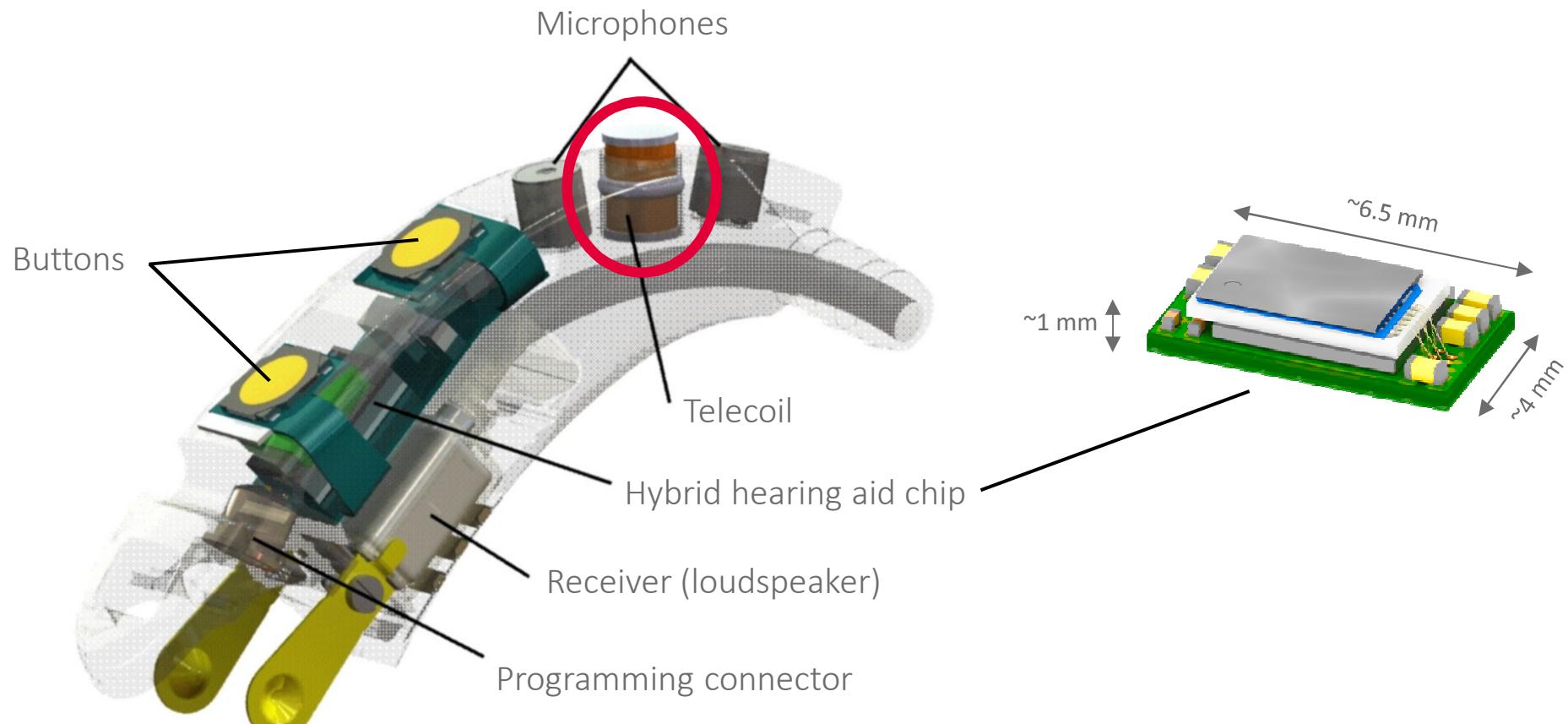

Source: kind.com, 2021


•• Excursion: Telecoil

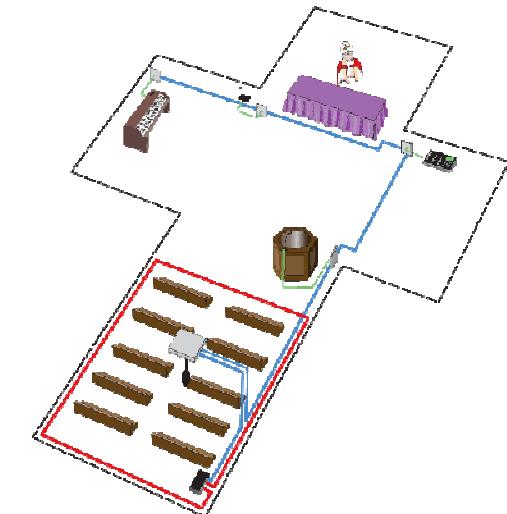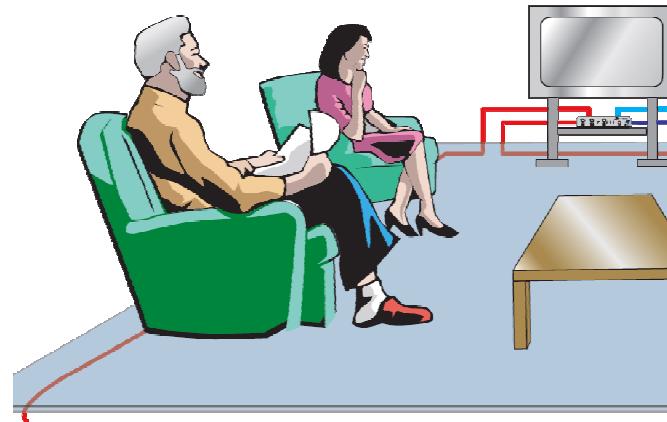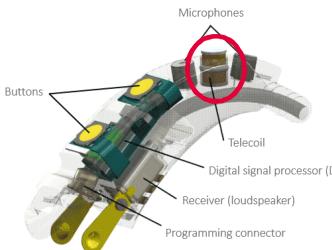



Source: signet-ac.co.uk, 2021

T-coil or induction loop receiver:


- Operating at audio frequencies, no RF.
- Induction loop parts:
 - constant current amplifier,
 - closed loop cable.
- Attention (privacy):
 - EM-field bleeding.
 - Always broadcast.

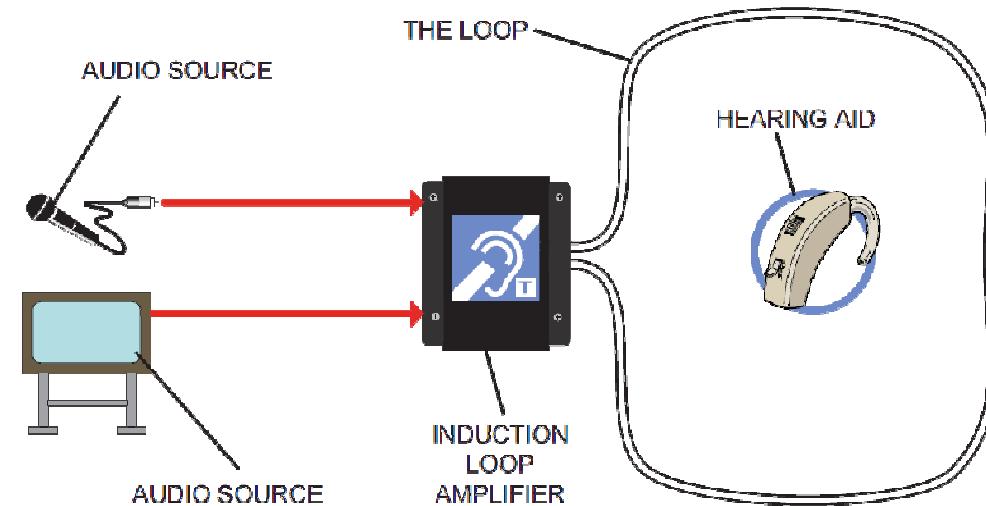

•• Recap: Equal-loudness contours

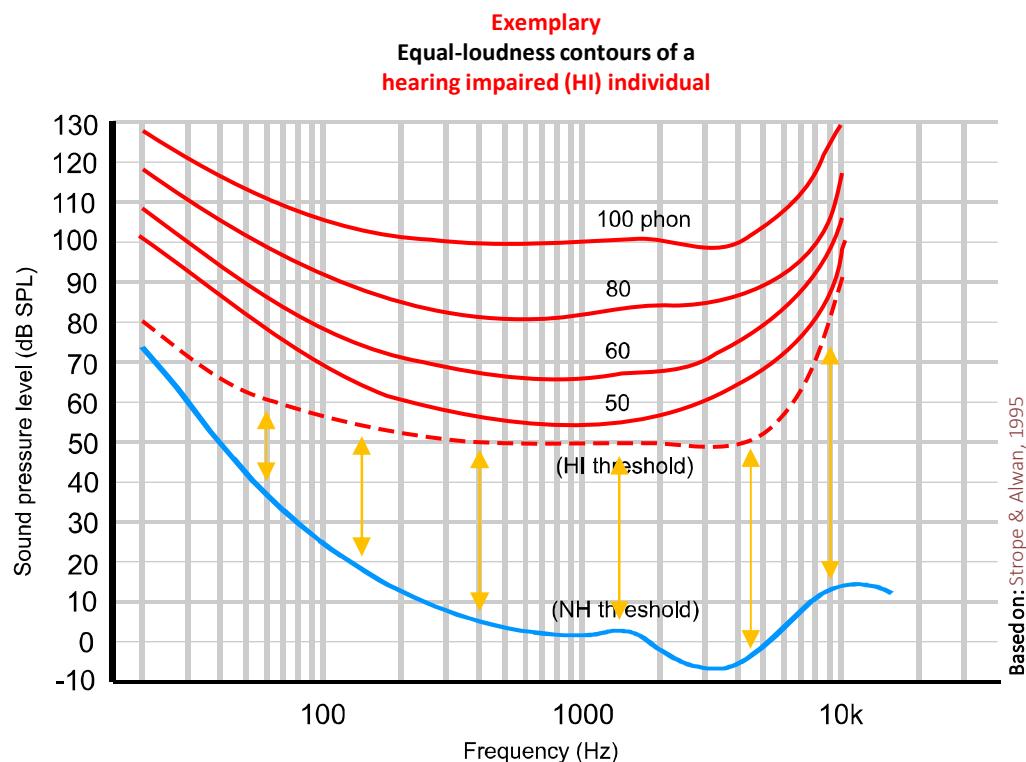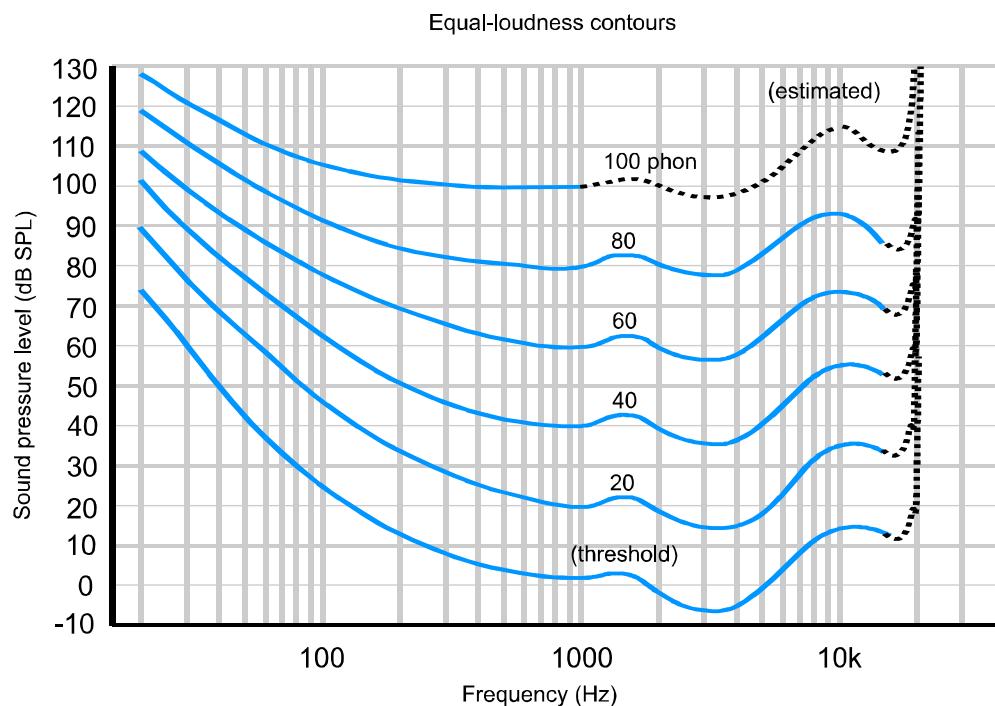


•• Hearing aids – Wide dynamic range compression

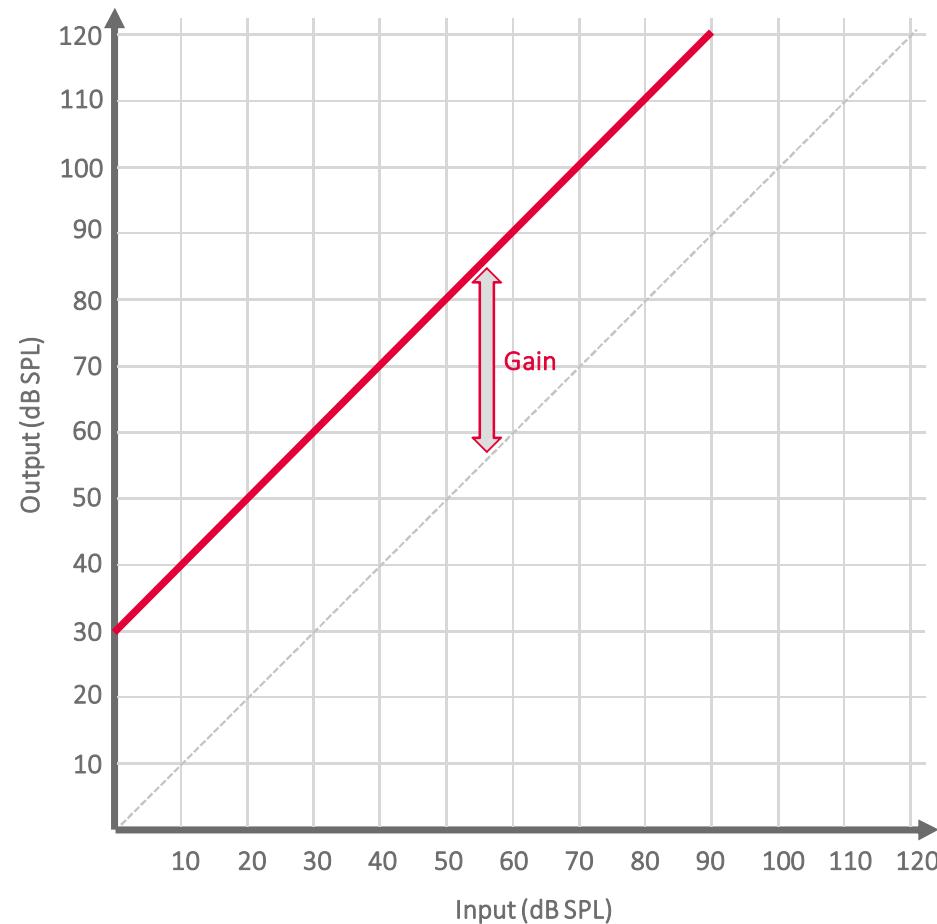

Linear amplification alone:

- Increases amplitude always with the same gain.
- Not conform with loudness perception of hearing impaired. (Will often cause uncomfortably loud percept.)
- Can lead to more hearing impairment.


•• Hearing aids – Wide dynamic range compression

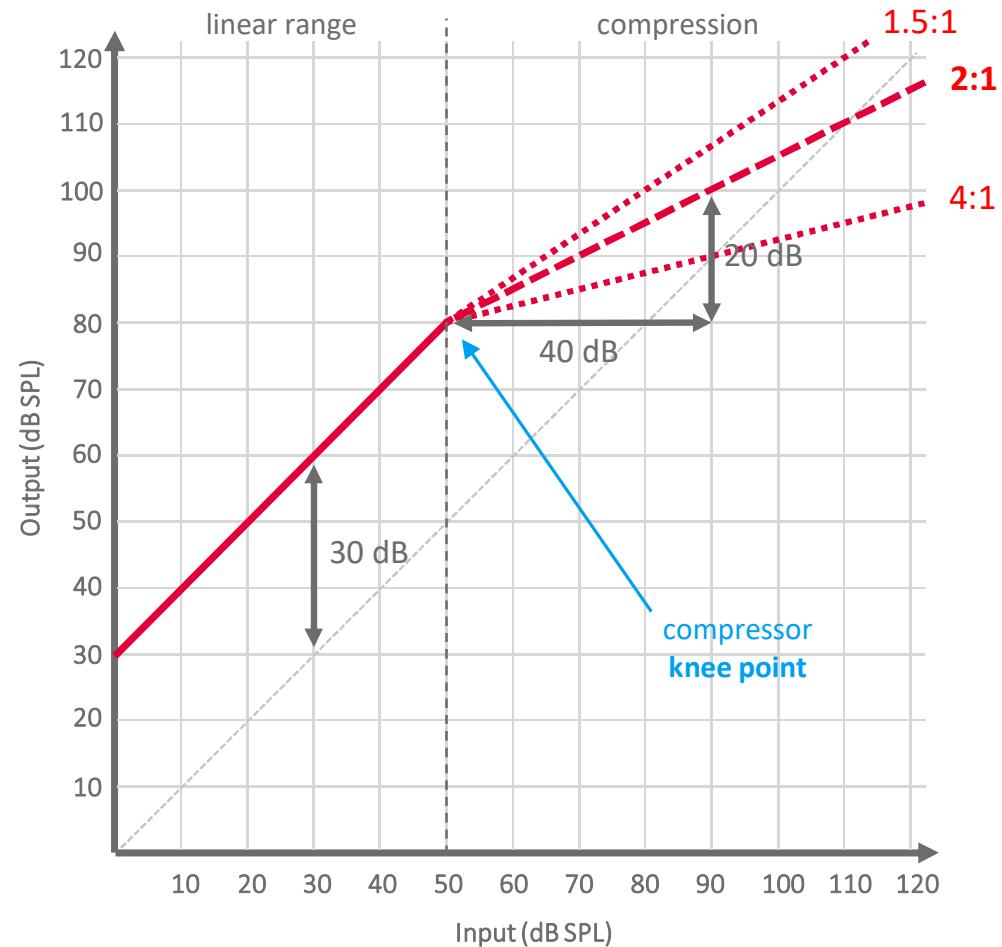

Compression:

- Example:
40 dB increase in the input leads to
20 dB increase in the output, i.e.,
→ compression ratio: $40 \text{ dB} : 20 \text{ dB} = 2:1$
- Compression rate can vary, but usually does not exceed 4:1 in the hearing aid field.


•• Hearing aids – Wide dynamic range compression

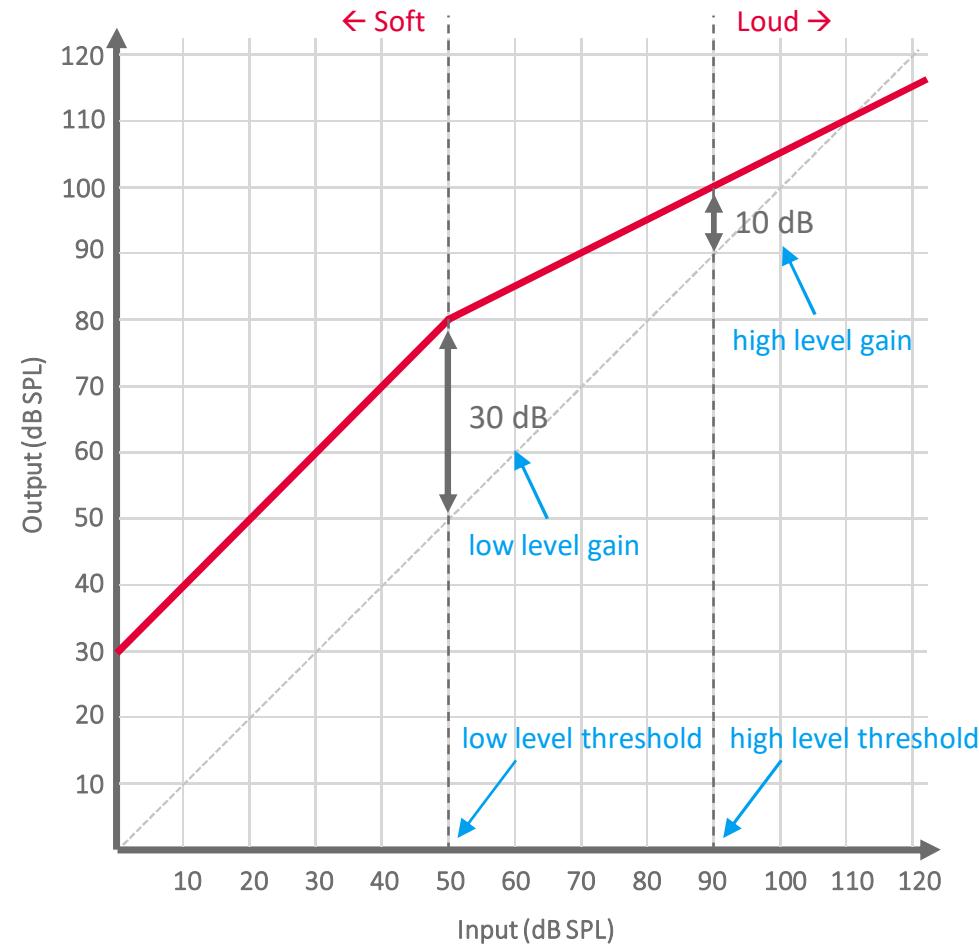

Compression:

- Low level gain: gain at ca. 50 dB SPL input
- High level gain: gain at ca. 90 dB SPL input


•• Hearing aids – Wide dynamic range compression

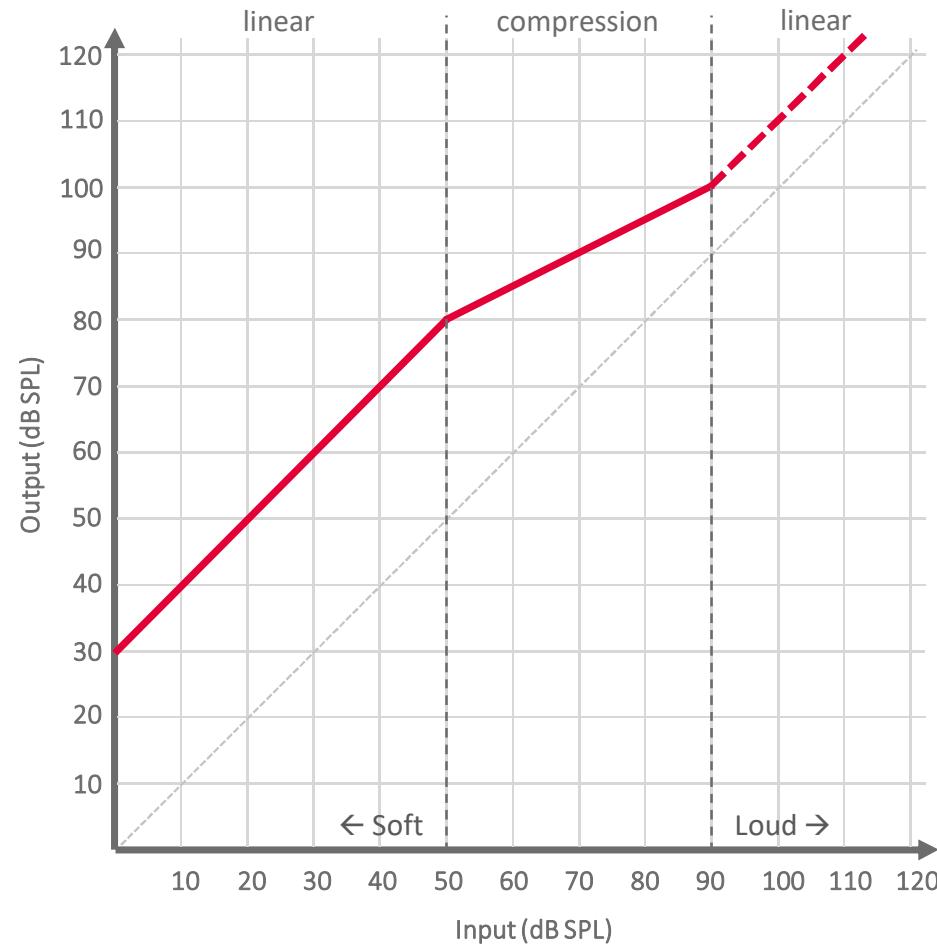

Second linear range:


- Can be useful to allow the listener to realize that the input is very loud.
- Should be verified to not extend into the threshold of pain. (See “limiter” on next slide.)
- Remember the normal hearing ear handles loud noises very similarly! Basilar membrane nonlinearity:

Recap

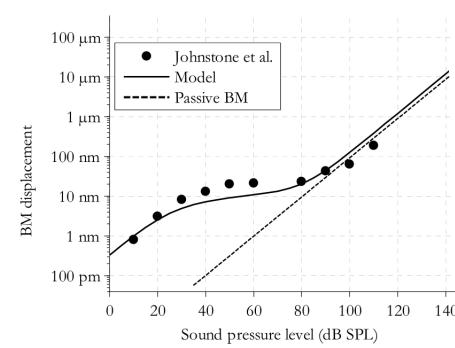

Adapted from: Johnstone et al., 1986


•• Hearing aids – Wide dynamic range compression

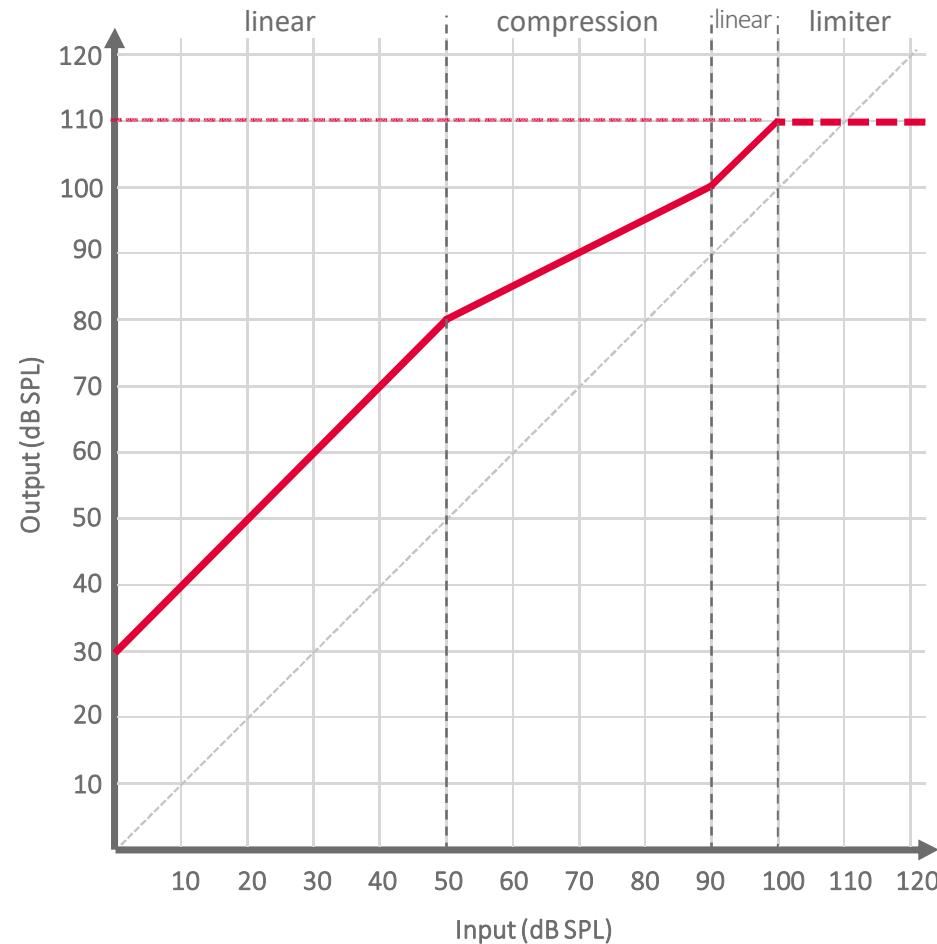

Limiter:

- Ensures not to cross the threshold of pain.
- Ensures technical limits not to be exceeded.
(MPO: maximum power output of a hearing device. Limited for example by the power source, amplifier, and receiver.)
- Is equivalent to a compressor with $\infty:1$ compression rate.


•• Hearing aids – Wide dynamic range compression

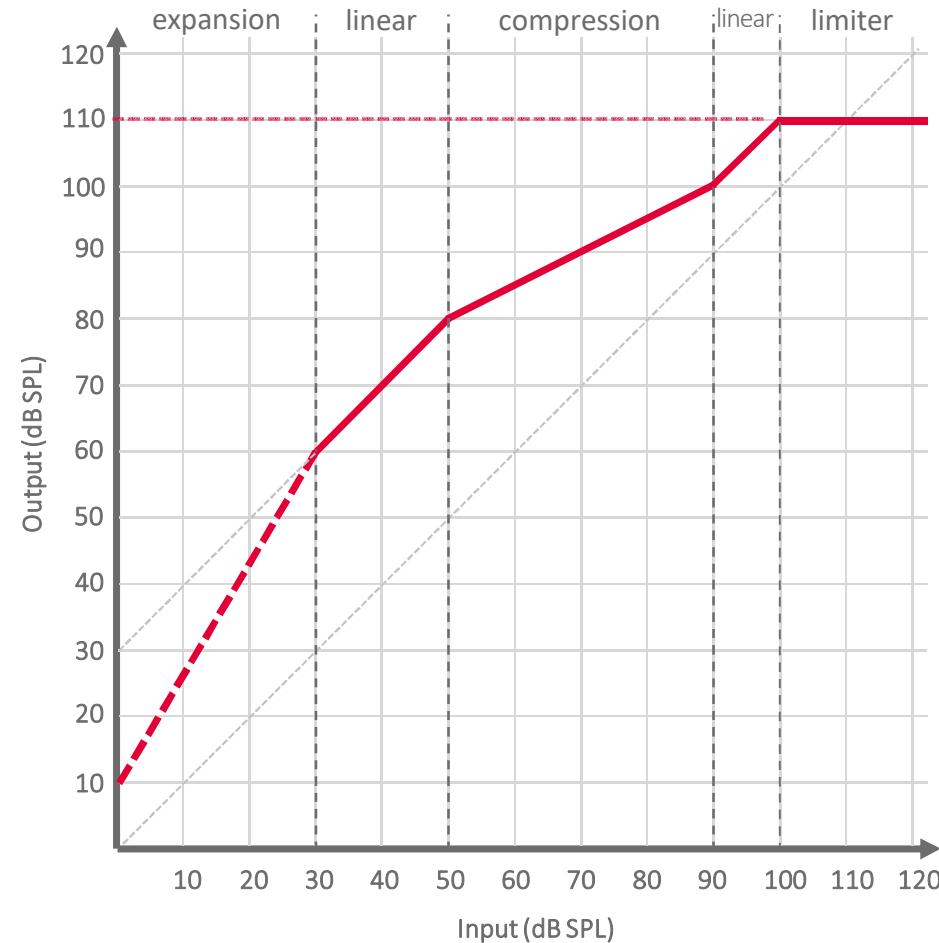

Expansion:


- Usually employed in the low SPL range.
→ Less gain added to soft noises.
- Reduces inherent noise of the system (e.g. that of microphones or preamps).
- Reduces irrelevant background noise.
→ Increases overall SNR.
- Typical expansion rate: 1:1.5 - 1:3.
(1:1.66 in the graph on the left.)

•• Hearing aids – Wide dynamic range compression

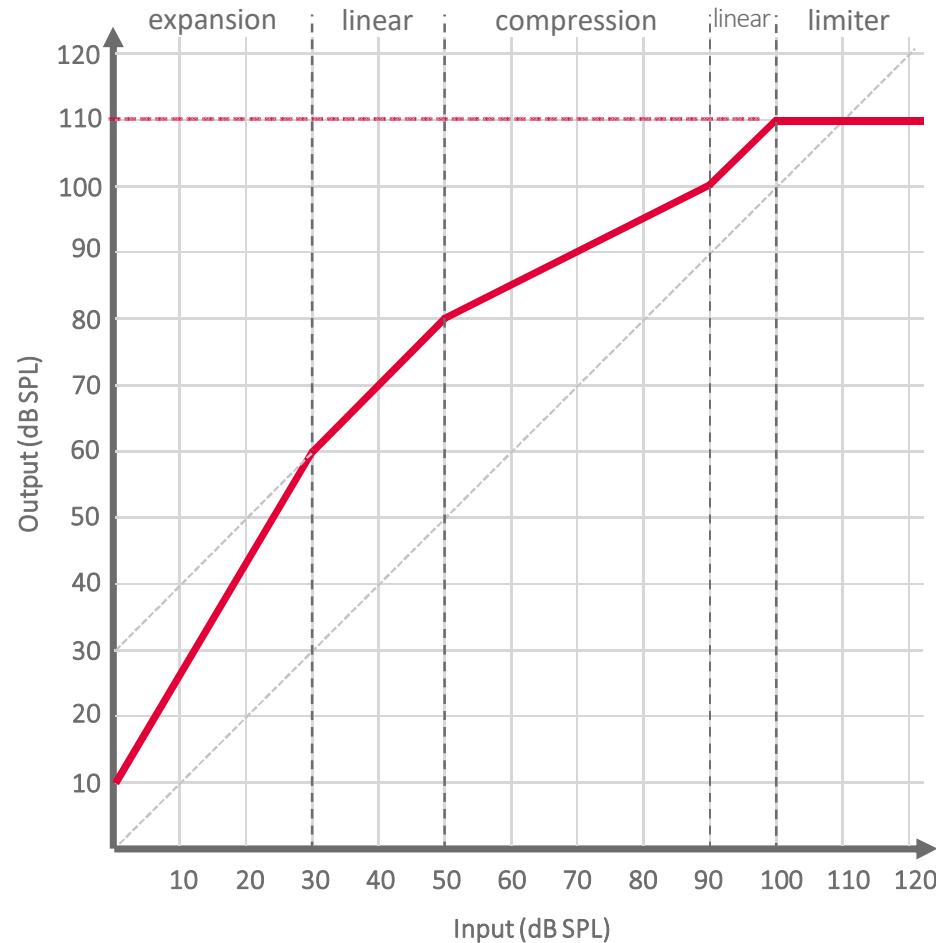


Gain plot:

- A shorter way to describe compression-expansion curves.

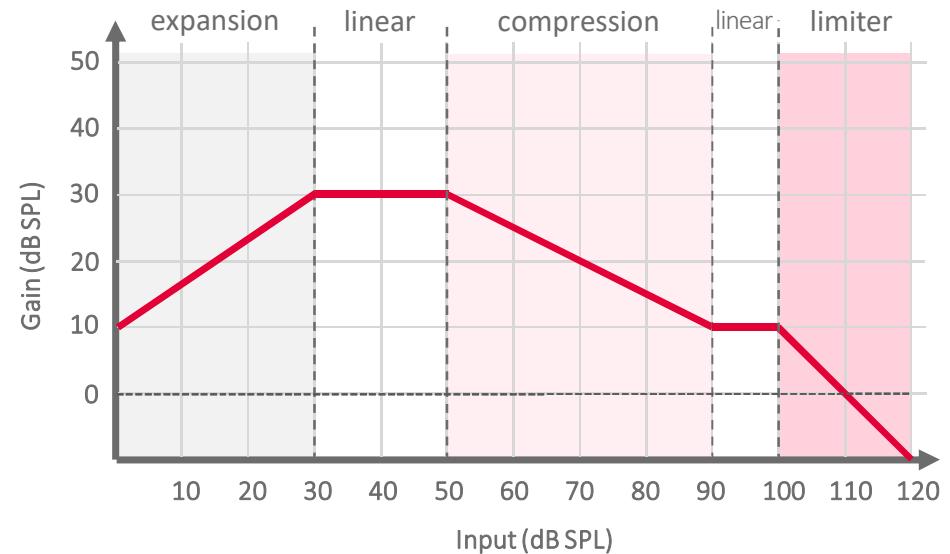


•• Hearing aids – Multi-channel WDRC

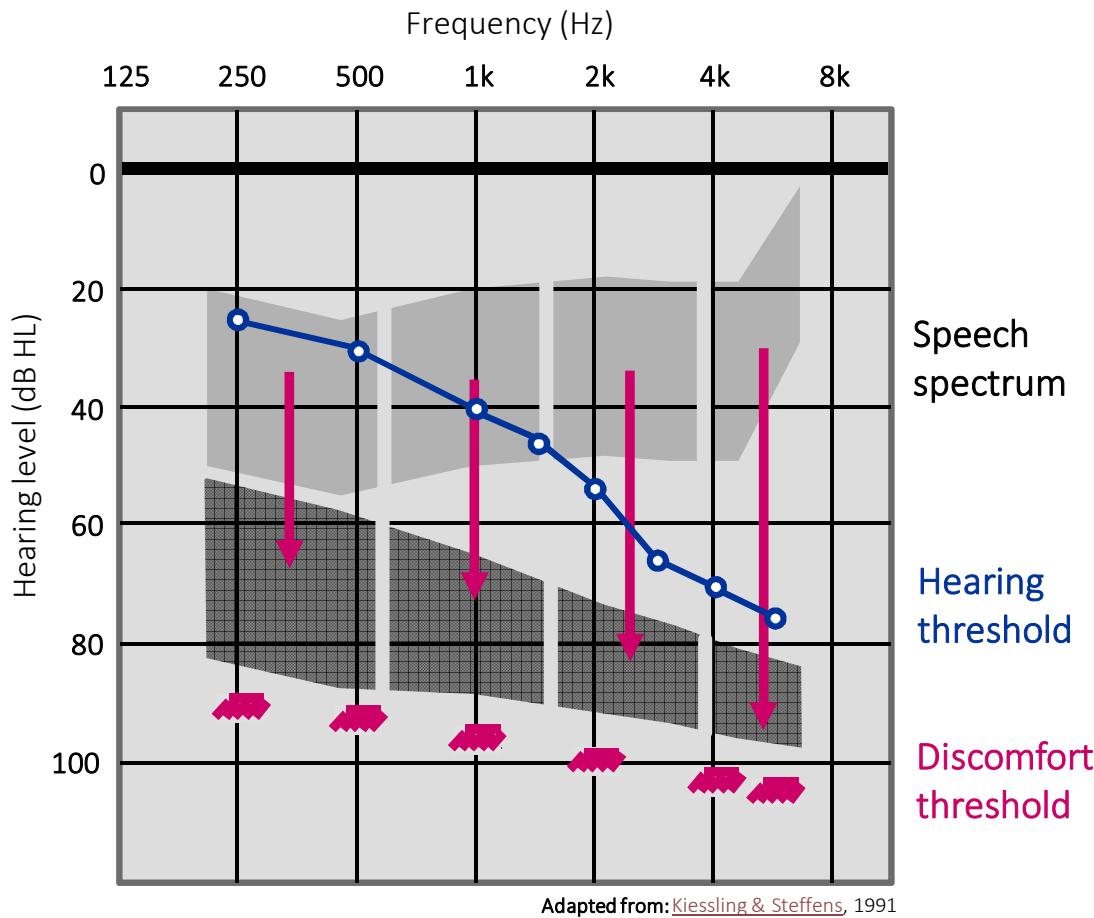


Multi-channel compression:

- The normal dynamic range of speech varies across frequency.
- The amount of hearing loss also varies across frequency, individually.
- Solution: use multiple dynamic compressors, each fitted for a distinctive frequency range.


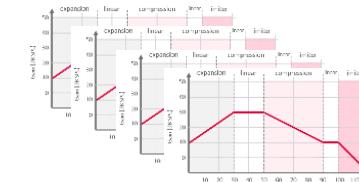


- A computationally cheap, but high quality filter-bank with almost perfect reconstruction is needed.
- (Keywords for searching further information: FFT-based WOLA filter-banks and warped filters.)


•• Recap: Causes of hearing impairment

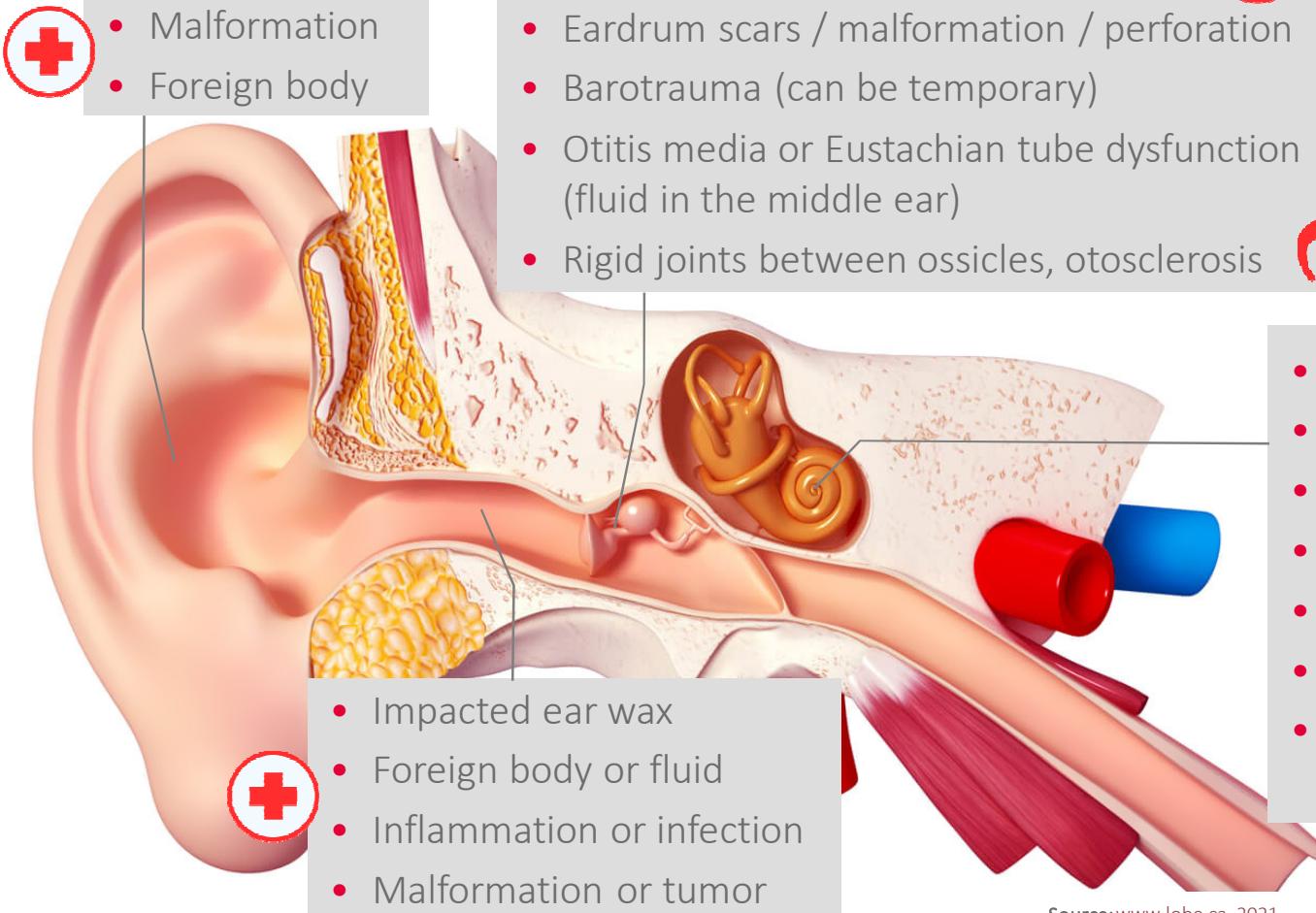

Source: www.lobe.ca, 2021.

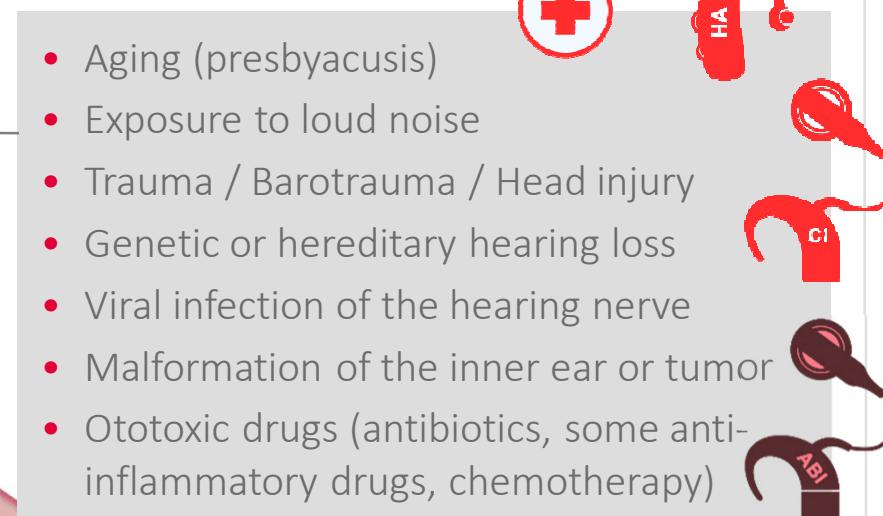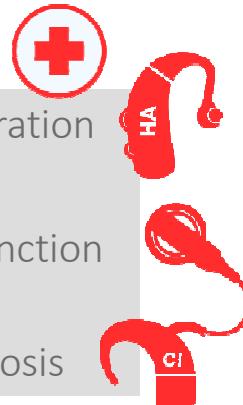

•• Cochlear implants

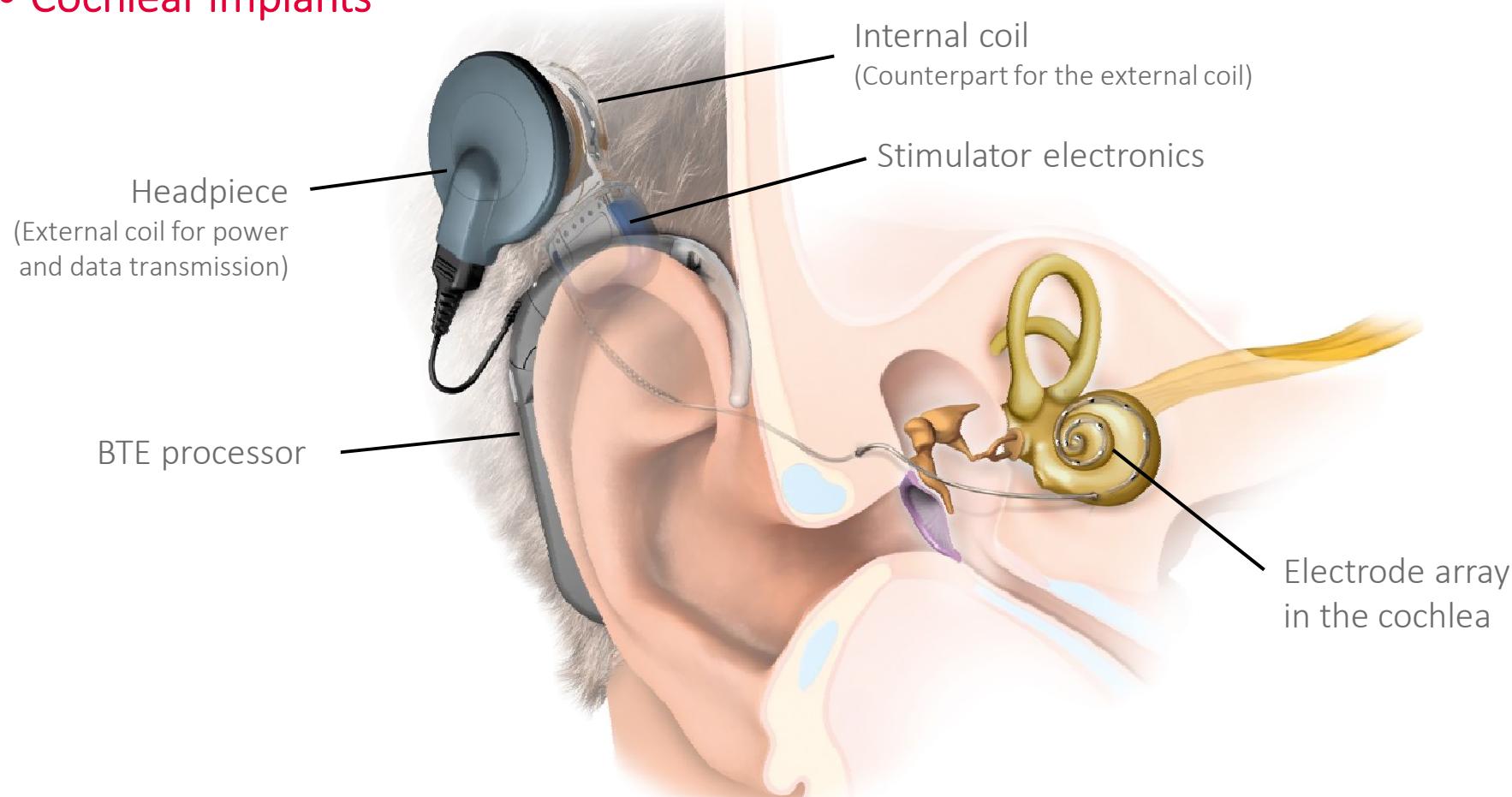

Source: medel.com, 2016

•• Cochlear implants




Source: [cochlear.com](https://www.cochlear.com), 2021 (audio changed via TigerCIS software for educational purposes only)

•• Cochlear implants



Source: cochlear.com, 2021 (unprocessed sound)


•• Cochlear implants – Signal processing strategy

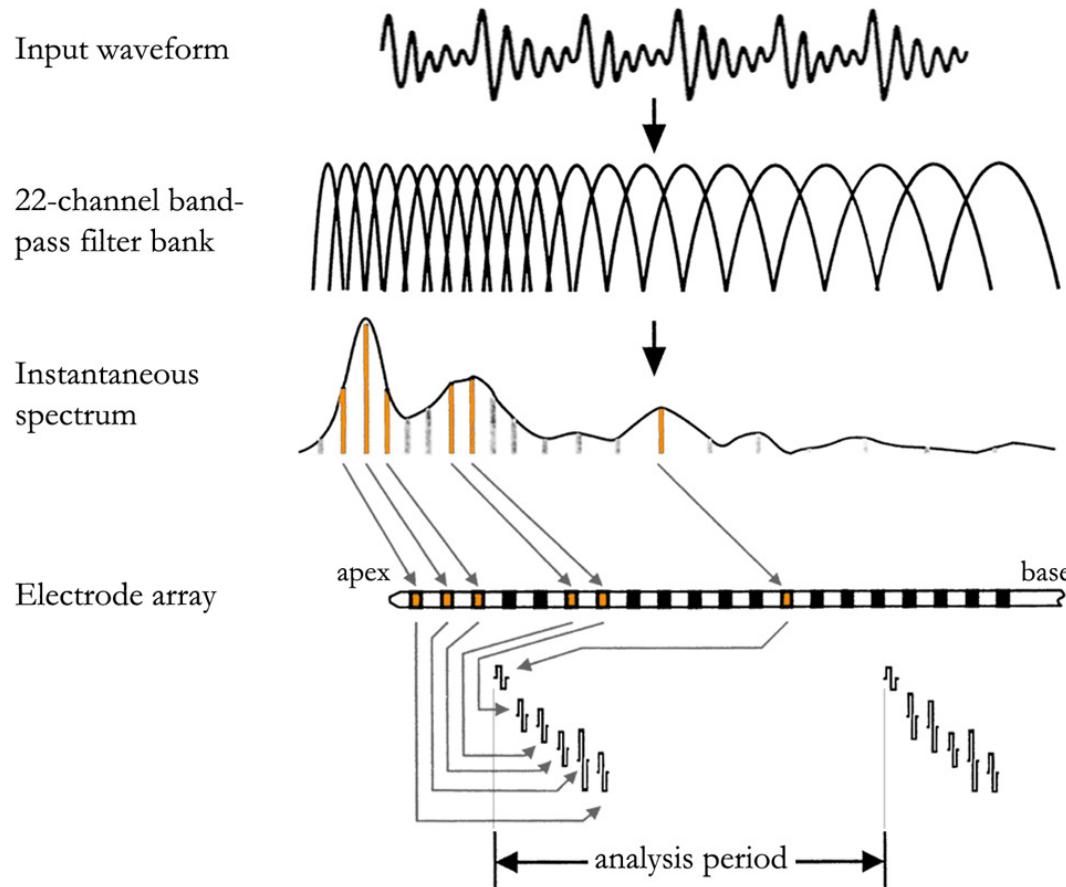

Source: [G. Clark, Cochlear Implants](#), 2003

ACE signal processing strategy
(advanced combination encoder)

Biphasic, charge balanced pulses

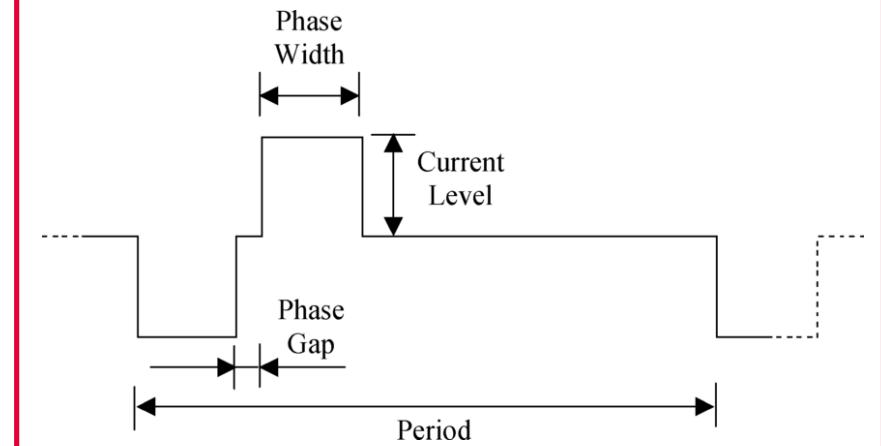

Source: [Tognola et al., 2005](#)

•• Cochlear implants – Early history

Source: hearinghealthmatters.org, 2021

Alessandro Volta
50V DC ear-to-ear



1790s

Source: abebooks.com, 2021.



Rudolf Brenner's
AC stimulation

Lundberg (1950),
Djourno & Eyries (1957)

Direct electrical stimulation of the auditory nerve with distinguishable percepts in human.


Source: Eisenberg.com, 2005



William F. House's
single-channel CI

House's first clinical
multi-channel CI

Source: Clark.com, 2013

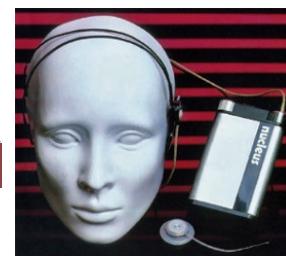

1961

1982

1985

1996

2001

Clarion implant
Advanced Bionics
(USA)



Source: cochlear.com, 2021

Nucleus implant
Cochlear Ltd. (AU)

Commercialization →

Combi40+ implant
MED-EL GmbH (AT)

•• Hearing aids (HA) vs. Cochlear implants (CI)

- System size

BTE processor is about the same size for both, but CI has more system parts (RF link, implanted stimulator, implanted electrode array).

- Power consumption

Signal processing requires about the same power for both, but CIs need extra current for electrical stimuli and RF link.

- Core signal processing

Both do spectral decomposition via an analysis filter-bank, but while HAs need to synthesize audio on the output, CIs calculate parameter for electrical stimuli.

- Additional algorithms / features

Very similar: directionality, dynamic range compression, noise reduction, binaural algorithms, but CIs do not need feedback canceller. Both can be streamed wirelessly and be remote controlled.

- Current issues

Both: speech in noise enhancement, power consumption and size, connectivity, user-friendly handling.
CIs: reverb, pitch discrimination (12/16/22/24 channels, does not matter ← current spread).

•• Cochlear implants – Ways to improve pitch sensitivity

- Train cochlear implant users

Use the brain's plasticity. Needed anyway in CI therapy, but has its limits.

- Enhance pitch changes


Smart preprocessing can increase speech melody changes in real-time (PREX project, Fraunhofer IDMT).



- Improve temporal pitch cues

Better stimulation strategies (SAM project, Fraunhofer IDMT)

•• Pitch perception

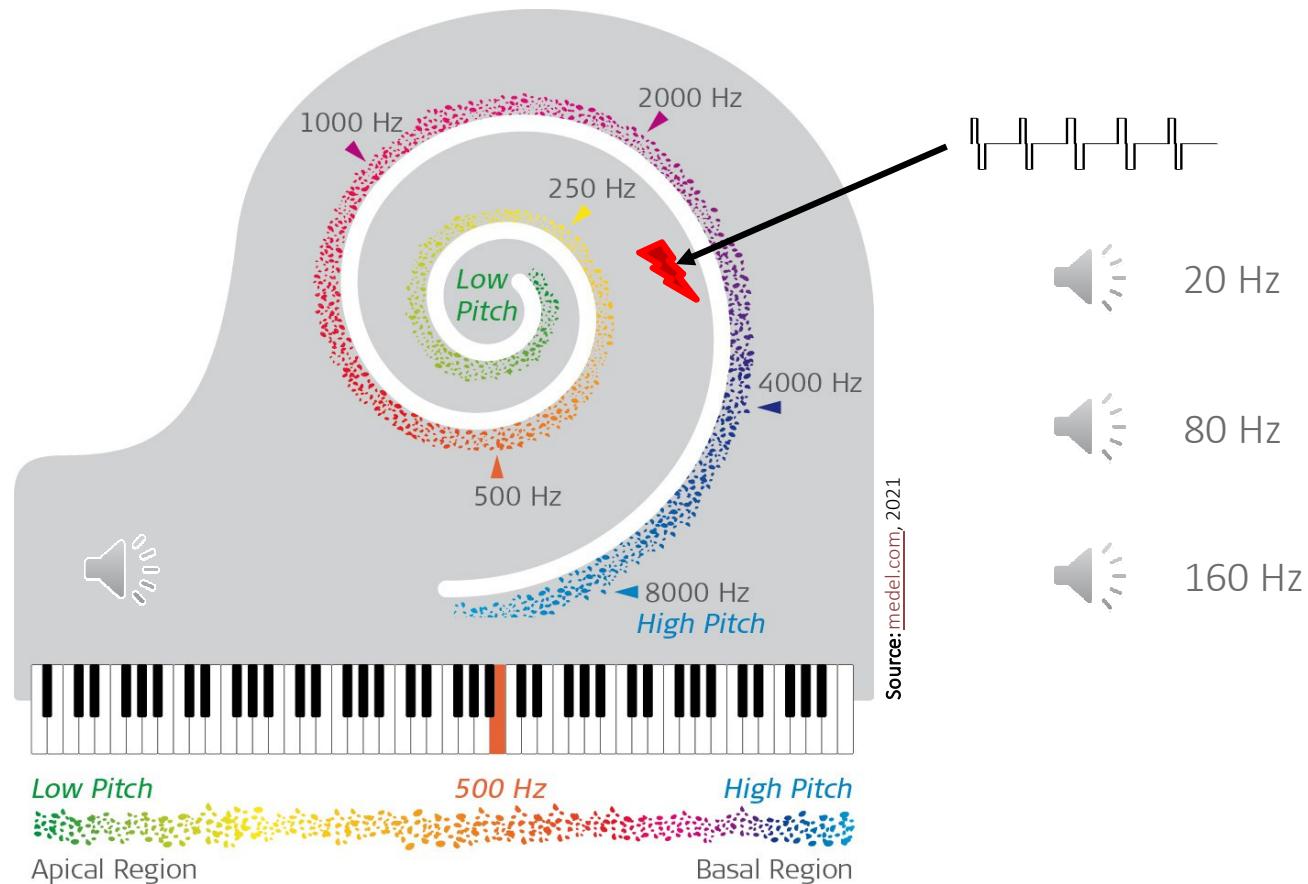

•• Cochlear implants – Ways to improve pitch sensitivity

- Train cochlear implant users

Use the brain's plasticity. Needed anyway in CI therapy, but has its limits.

- Enhance pitch changes

Smart preprocessing can increase speech melody changes in real-time (PREX project, Fraunhofer IDMT).



- Improve temporal pitch cues

Better stimulation strategies (SAM project, Fraunhofer IDMT)




- Improve place pitch cues (reducing current spread)


- Better electrode design ("modiolus-hugging" or "adhering" electrodes),
 - Current steering (parallel stimulation over multiple electrode contact points)
 - Optogenetic stimulation (optical cochlear implants).

•• Optogenetic cochlear implants


Electrical stimulation

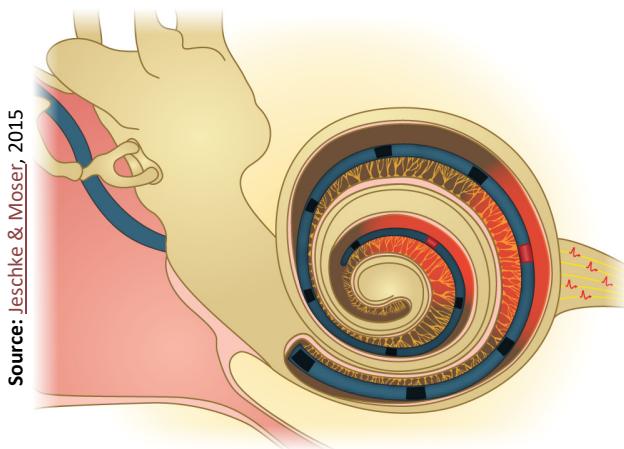

Optical stimulation

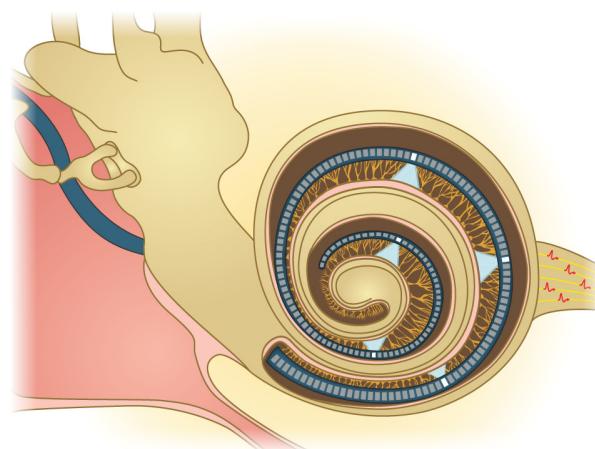

Research hardware

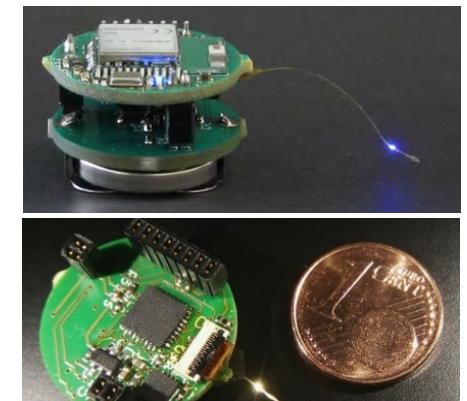


Challenges of optogenetic cochlea stimulation:

- **Light-sensitive neurons:** efficient, reliable and safe opsin* expression via *genetic manipulation* of the spiral ganglion neurons.
- **Stimulation technology** must be efficient, reliable and safe.


*opsin = light-gated ion channel

•• Optogenetics

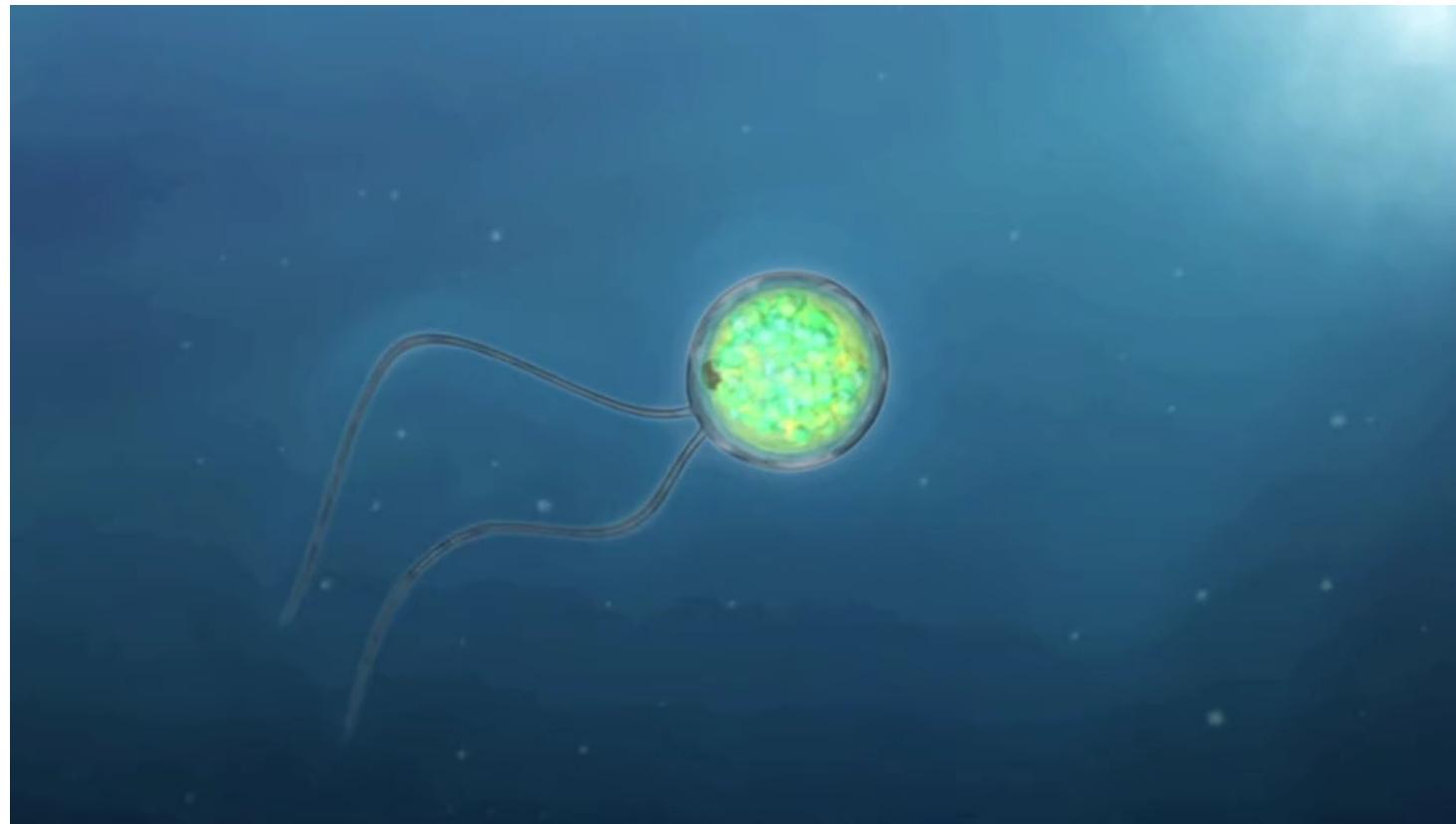

Source: be.mit.edu, 2021

- Optogenetics – Demo with the research hardware

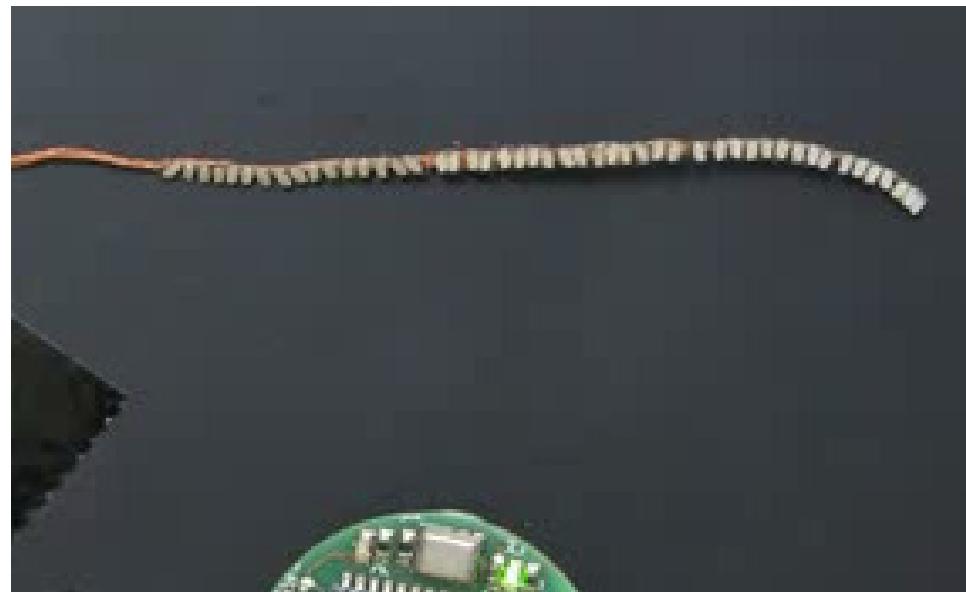

Music: [Josh Armistead - Lazer Beams](#)

Licence: Attribution 4.0 International (CC BY 4.0)



Thank you very much!
Questions?

Dr.-Ing. Tamas Harczos
tamas.harczos@audifon.com

audifon GmbH & Co. KG